

# CiliaGUI – Installation and User Guide

---

## Installation

This video describes the installation of the CiliaDetect App.

To run the toolbox, you need to:

1. Install MATLAB (2021 or later).

- MATLAB is not free; check with your institution for access.
- If you cannot get a license, contact me — we may try a standalone version of the program.
- Installation instructions: <https://www.mathworks.com/help/install/ug/install-products-with-internet-connection.html>
- During installation, make sure to install the following toolboxes:

- **Statistics and Machine Learning Toolbox**

- **Image Processing Toolbox**

2. Download the code from GitHub:

- Repository: <https://github.com/pmgarderes/CiliaDetect2.0>
- Click 'Download as ZIP file'.

3. Unzip and place the folder on your hard drive.


- Example: on Windows, I use WinRAR and put the files into a folder named ciliaGUI\_demo at the root of C:/.

## Launching a Session

4. Navigate to the folder and find Full\_GUI.m.

- Double-click it. If MATLAB is not open, this will automatically launch MATLAB and set the current folder correctly.

- Make sure the current folder in MATLAB (shown in the top bar) contains Full\_GUI.m. If it doesn't, adjust it manually. → it should look like this (roughly)

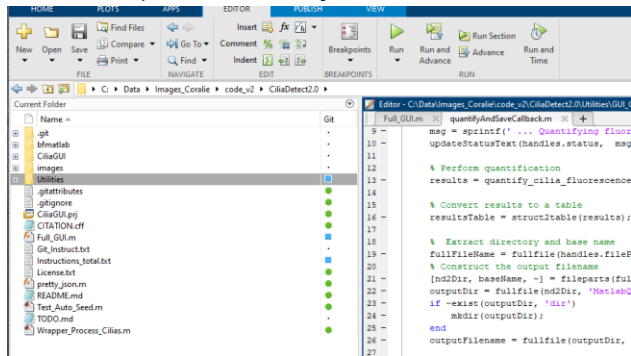

## 5. Start a session:

- Go to the Editor panel.
- With Full\_GUI.m active (double click Full\_GUI in matlab to make sure it's active)
- click Run.

## Preprocessing – Z-Planes

### 6. Prepare your microscopy images before using the GUI.

- The program only reads \*\_reduced.mat, i.e. files formatted for matlab. Thus we need to pre-process microscopy image and format them.
- Strategy: reduce the number of frames in a to a few averaged Z-planes (down to a single Z-plane), so multiple Z-planes can be separated from the same Z-stack.

To preprocess:

1. Click the 'Pre-process' button.
2. Enter the downsampling factor (number of frames to average together).
  - Example: enter '1000' to collapse the full stack into a single Z-plane.
3. A file browser window will open. Navigate to the folder containing your microscopy files.
  - You will not see the individual files — just select the directory where the files are located.
4. The GUI text box will display messages ('WAIT') and begin downsampling.

Notes:

- Preprocessing large stacks may take a long time. Tip: you can batch-process many files overnight.

- Output: a new folder 'Matlab\_quantif' is created. Each original file generates a reduced version with the suffix '\_reduced.mat'.

## Loading a File

7. Load files into the GUI:

- Only files ending in '\_reduced.mat' can be loaded.
- Click 'Select File', navigate to the directory, and choose the file.

## GUI Display

8. The central panel displays one Z-plane in one channel.

- Adjust fluorescence intensity with numeric keypad: /, \*, -, +.
  - These settings apply across all Z-planes in the same channel but not across channels.
- Navigate between channels: left/right arrow keys.
- Navigate between Z-planes: up/down arrow keys.

## Detecting ROIs

9. To add a new ROI:

- Left-click directly on the cilium (not beside it).
- Press Space.
- The program segments the cilium around your click.
- If detection fails (e.g., low contrast), it will appear as a single point.
  - Delete failed detections by pressing U once (undo).

☐ Do not leave empty ROIs.

## Parameters to Adjust Detection

10. Sensitivity is the key parameter.

- Recommended range: 0.3–0.7 (use 0.4–0.5 in most cases).

- Higher sensitivity: easier detection, less accurate ROIs.
- Lower sensitivity: stricter detection, fewer but more accurate ROIs.

## Merging ROIs

11. If a cilium is split across Z-planes or incompletely detected:

- Click near the junction point.
- Press M to merge ROIs (works if they are close enough).

## Deleting ROIs

12. To delete ROIs:

- Single ROI: click near it and press S.
- All ROIs: click 'Clear Detections'.

## Other Useful Buttons

13. Refresh ROI display: press R.

## Saving Detections

14. Save your work:

- Click 'Save Detections' at any time.
- This creates a new file in 'Matlab\_quantif' inside 'reduced\_stack/'.
- File name format: originalName\_reduced\_cilia\_detections.mat
- Next time you open the reduced file, saved detections load automatically.

☐ Only load '\_reduced.mat' files. Do not load detection files directly.

## Visualizing Masks

15. To view ROI and background masks:

- Click 'Visualize Masks'.
- Two images will appear:
  - Left: cilia mask overlay.

- Right: background mask.
- Images are saved automatically.
- Adjustable parameters:
  - Padding: gap between ROI and background (pixels).
  - Background mask width (pixels).

## Quantify and Save Fluorescence

16. To export quantified data:

- Click 'Quantify and Save Fluorescence'.
- Creates one Excel file (.xls) per image in 'Matlab\_quantif/'.
- Each row = one cilium.

Fluorescence outputs (per channel):

- F\_StackSum\_ch# – sum across the Z-stack.
- F\_plan\_ch# – mean in the Z-plane.

Morphological measures (if metadata available):

- Length ( $\mu\text{m}$ )
- Width ( $\mu\text{m}$ )
- Curviness (length/chord length)
- Length/Width ratio
- Area ( $\mu\text{m}^2$ )

Notes:

- Both methods correct for background.
- Both are valid for comparisons if exposure and Z-range are similar.
- Recommendation:
  - Use stack sum if Z-range > cilium span.
  - Use mean method otherwise.

## Editing Parameters

17. Parameters are user-editable:

- Any parameter can be modified.
- Settings are automatically saved and reloaded in the next session.